


SETTING UP YOUR KNEE REST SUPPORT BRACE

The support brace is intended to be used to retain your Astor Bannerman knee support in place; it must NOT be used for any other purpose.

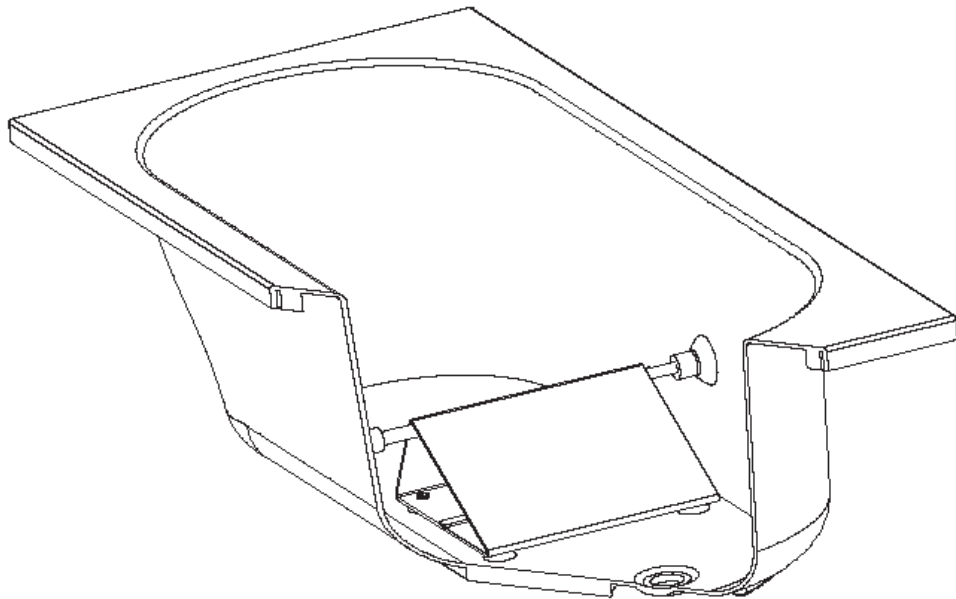

TELESCOPIC BRACES

The Support brace is telescopic – it can easily be adjusted to fit most baths

Slide the support through the knee support.

Lubricate the suction cups on the end of the back brace with mild soap and water

Extend the back brace to fit the sides of the bath, ensuring they are pushed firmly against the sides of the bath.

Please check that the suction cups are firmly in position before using the knee support.

SUPPORT REMOVAL

To detach the support from the bath, lift an edge of each suction cup to release the pressure and remove the telescopic support.

Lift the edge of each suction cup to break the seal and remove the knee support.

CARE AND MAINTENANCE

N.B Cleaning is recommended on a regular basis

The brace & suction cups should be cleaned after use with mild soap and water.

Rinse with clean water, and dry.

Store the chair in a cool dry place out of direct sunlight.

DO NOT use bleach, solvents, abrasives, synthetic detergents, wax polishes, antibacterial sprays or wipes.

Always keep this product away from naked flames, cigarettes and sources of heat including open fires, radiators and heaters.

For further information please refer to MHRA or your NHS cleaning guidelines.

DAILY CHECKS

Check the product for loose, damaged or missing components

Check the brace can be adjusted correctly, and that the suction cups firmly locate the knee support in the bath

Keep all parts clean.